

DtM-v4 Headband Strap


8 April 2020

Tests:


This variety of strap used to secure the face shield HEADBAND FRAME to the care provider's head was tested using four common cleaning solutions. This headband strap was made of non-latex silicone rubber MSDS CW9003. Soap and water, water with 10% bleach, Super Sani-Wipes, and Cavi-Wipes were separately wiped on to the surface of the bands. The water with soap and the water with 10% bleach tests were wiped dry. The Super Sani-Wipes and the Cavi-Wipes were left to air dry. This process was repeated for each cleaning solution five times.

Results:

The surface of the HEADBAND STRAP was easy to clean. There was no discoloration or change in the strength properties of the band material from cleaning.

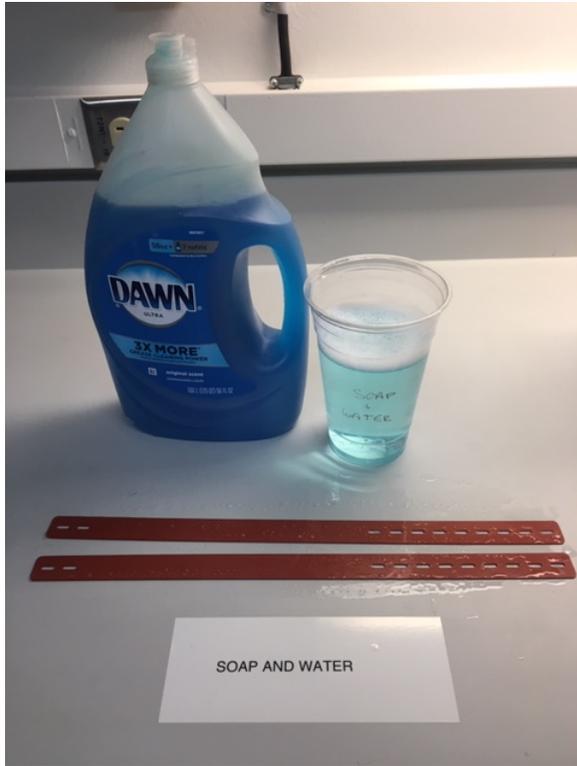

SOAP AND WATER

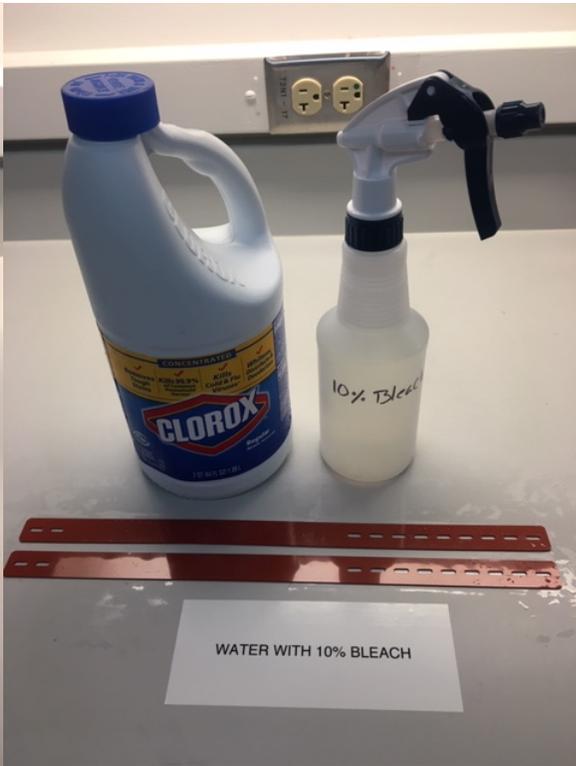

WATER WITH 10% BLEACH



SANI-WIPES



CAVI-WIPES